

MJB - Under Bed Light (ext.), Out of Bed Detection



To learn more about LINAK please visit:
WWW.LINAK.COM

LINAK® 
WE IMPROVE YOUR LIFE

Contents

Preface	4
Important information	5
Safety instructions	5
Residual risk	6
Before installation, reinstallation or troubleshooting	7
Before start-up	7
During operation	7
Information on the MJB	7
Repairs	8
Misc. on the MJB system	9
Warranty and service life	9
Maintenance	9
MJB with UBL (ext.)	10
Standard functionality	10
Features	11
Options (speciale article)	11
Usage	11
Precautions	12
IP Protection degree	13
Environmental conditions	14
Modular plug with wide and narrow slit	15
Keys to symbols	16
Drawing appendix	18
Addresses	20

Preface

We are delighted that you have chosen a product from LINAK®.

LINAK systems are high-tech products based on many years of experience in the manufacture and development of actuators, electric control boxes, controls and chargers.

This User Manual will tell you how to install, use and maintain your LINAK equipment.

We are sure that your LINAK system will give you many years of problem-free operation. Before our products leave the factory they undergo full function and quality testing. Should you nevertheless experience problems with your systems, you are always welcome to contact your local dealer.

LINAK subsidiaries and some distributors situated all over the world have authorised service centres, which are always ready to help you.

LINAK provides a warranty on all its products. This warranty, however, is subject to correct use in accordance with the specifications, maintenance being done correctly and any repairs being carried out at a service centre, which is authorised to repair LINAK products.

Changes in installation and use of LINAK systems can affect their operation and durability. The products are not to be opened by unauthorised personnel. The User Manual has been written on the basis of our present technical knowledge. We are constantly working on updating the information and we therefore reserve the right to carry out technical modifications.

LINAK A/S

Important information

Description of the various signs used in this manual.

Warning!



Failure to comply with these instructions may result in accidents involving serious personal injury.



Failing to follow these instructions can result in the product being damaged or destroyed.

Safety instructions



Safe use of the system is possible only when the operating instructions are read completely and the instructions contained are strictly observed.

Failure to comply with instructions marked with the "NOTE" symbol may result in serious damage to the system or one of its components.

Persons who do not have the necessary experience or knowledge of the product/products must not use the product/products. Besides, persons with reduced physical or mental abilities must not use the product/products, unless they are under surveillance or they have been thoroughly instructed in the use of the apparatus by a person who is responsible for the safety of these persons.

Moreover, children must be under surveillance to ensure that they do not play with the product.



It is important for everyone who is to connect, install, or use the systems to have the necessary information and access to this User Manual.



If there is visible damage to the product it must not be installed.



The appliance is not intended for use by young children or infirm persons without supervision.



Young children should be supervised to ensure that they do not play with the appliance.

**Residual risk**

Some of the products contains software based components. LINAK has made various possible efforts to assure that the software is free of errors and that the software has been developed according to the rules of IEC 60601-1-4 (software in Medical products). That involves Risk Analysis which shows a small residual risk for unwanted/unintended movement of actuators under specific conditions.

According to the above rules it must be informed and if necessary considered in the Risk analysis of the final application - More details to residual risk can be provided by LINAK if necessary.

**Warning!**

The plastic parts in the system cannot tolerate cutting oil.

**Warning!**

LINAK's actuators and electronics are not constructed for use within the following fields:

- Planes and other aircrafts
- Explosive environments
- Nuclear power generation

**Warning!**

If faults are observed, the products must be replaced.

**Warning!**

Never spray directly on the products with a high-pressure cleaner.

**Warning!**

Always insure that the patient is accessible in case of emergency.

**Warning!**

Cables must be secured safely so that persons cannot get trapped in it.

**Warning!**

Always use locking mechanism and O-ring

**Warning!**

Sockets not used must be fitted with blind plugs to ensure IP Degree



The LINAK products cannot tolerate the influence of strong solvents, basic or alkaline liquids.

Before installation, re-installation, or troubleshooting

- Stop the application
- Switch off the power supply and pull out the mains plug.
- Relieve the application of any loads, which may be released during the work.

Before start-up

- Make sure that the system has been installed as instructed in the User Manual.
- Make sure that the voltage of the control box is correct before the system is connected to the mains.
- System connection. The individual parts must be connected before the control box is connected to the mains.

During operation

- Ensure that the cables are not damaged.
- Unplug the mains cable on mobile equipment before it is moved.

Information on the MJB

The MJB Modular Junction Box is designed for use together with OpenBus™ control boxes.

The MJB makes it possible to connect multiple handsets, attendant controls or it makes it possible to use the MJB as a control unit for 3rd party products as Out of Bed detection, WET sensor, USB charger, Under Bed Light etc.



LINAK A/S delivers only the OpenBus system and is not responsible for any products other than LINAK products (i.e. products from 3rd party suppliers) or the compliance of such products with the LINAK OpenBus™ system.

Repairs

In order to avoid the risk of malfunction, all repairs must only be carried out by authorised LINAK workshops or repairers, as special tools must be used and special gaskets must be fitted. Products under warranty must also be returned to authorised LINAK workshops.



Warning!

If any of the LINAK products are opened, there will be a risk of subsequent malfunction.

Disposal of LINAK's products

LINAK's products may be disposed of, possibly by dividing them into different waste groups for recycling or combustion.

We recommend that our product is disassembled as much as possible at the disposal and that you try to recycle it. As an example of main groups within sorting of waste we can mention the following:

Metal, plastic, cable scrap, combustible material and collection for recoverable resources.

Some of these main groups can be sub-divided into groups e.g. metal can be divided into iron, stainless steel and aluminium and alloy steel. Plastic can e.g. be divided into ABS, PA, PE and PP.

The MJB can be disassembled and sorted into the following recycling groups:

Component	Recycling group
MJB	Plastic housing – Plastic recycling or combustion
	Cable – Cable scrap or combustion
	PC-board – Electronics scrap

By now all our casted plastic units are provided with an interior code for plastic types and fibre contents, if any.

Misc. on the MJB system

Warranty and service life

The LINAK warranty covers manufacturing defects in the products, starting from the date of manufacture. There is 18 months' warranty for MEDLINE and CARELINE products. The warranty is limited to the value of the LINAK product.

LINAK's guarantee is only valid so far as the products have been used and maintained correctly and has not been tampered with. Furthermore, the products must not be exposed to violent treatment. In the event of this, the warranty will be ineffective / invalid. LINAK's warranty is only valid if the system is unopened and has been used correctly.

All LINAK products are designed to have an optimum service life as a matter of course, but the expected service life in a specific application is very dependent on how the products are used.

Maintenance

- The LINAK products must be cleaned at regular intervals to remove dust and dirt and inspected for mechanical damage, wear and breaks.
- The LINAK MJB is a closed unit and requires no internal maintenance.
- IPX6 is waterproof.
- Make sure that the plugs in the MJB are correctly fitted with O-rings before washing.
- Electronics must be inspected at attachment points, wires, cabinet and plugs.
- Inspect the connections, cables, cabinets and plug, and check for correct functioning.


MJB with UBL (ext.), Switch input (S1/S2), version 506-01 & 506-02 :

The MJB is an MJB with 3 different options, Under Bed Light (UBL), switch input S1 which can be used for an external switch and switch input S2 to be used for Out of Bed Detection (OOBD).

UBL: The MJB with Under Bed light (UBL) is a simple solution to prevent falling accidents and make the patient feel safe. The UBL makes it easy for the patient to find their way back to the bed at night without disturbing the other patients. The MJB has an external LED cable (0964135) which is connected to Port 5. The external LED cable makes it more flexible to use the UBL. It can be moved from side to side of the bed or it can be placed at the foot end of the bed.

External Switch (S1): It is possible for the customer to connect a switch directly to the MJB. This could be used with a customised switch or control.

OOBD (S2): The Out of Bed Detection option is used when an alarm is needed when the patient leaves the bed. The OOBD is used in combination with a 3rd. party product, please refer to the 3rd party accessories area on the OpenBus™ site on LINTRA.

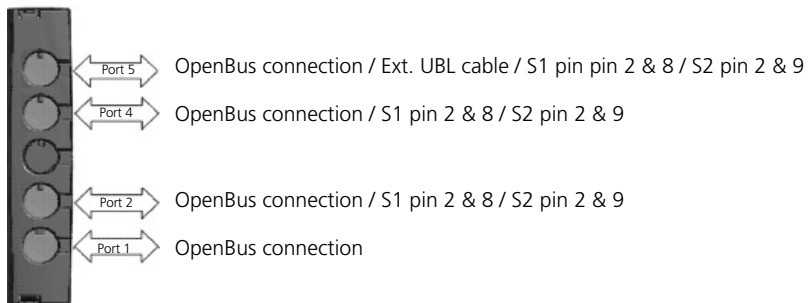

Standard functionality:

UBL: The UBL LED can be switched on and off via the patient control (Key1) or the attendant control (Key2).

External Switch (S1): The external switch is connected to Pin 2 & 8 on Port 2, Port 4 or Port 5.

The switch input functionality can be enabled/disabled via the attendant control. The enable/disable Status (switch status) is indicated on the OpenBus™.

The switch input functionality is as standard to be used with a NO switch. When the switch is being activated (NC), an alarm is being sent on the OpenBus™ (switch alarm).

Note! It is possible to combine OOBD, Switch input and OpenBus™ connection on the same port.

OOBD (S2): The OOBD sensor is connected to Pins 2 & 9 on Port 2, Port 4 or Port 5. The OOBD functionality can be enabled/disabled via the attendant control. The enable/disable status (OOBD status) is indicated on the OpenBus™. The OOBD function is as standard to be used with a NO (Normally Open) sensor (the sensor is NC (Normally Closed) when the patient is in the bed). When the sensor status is changing state from NC to NO (patient leaves the bed) the MJB will send an alarm (OOBD alarm) on the OpenBus™. The OOBD is to be used together with a Gateway MJB.

Features:

- Add-on value for hospital & care bed
- Compact and functional design which allows easy mounting
- 4 connections
- Connection to OpenBus™ control box via a modular plug cable, 6 wires
- Colour: light grey RAL 7035
- Protection class IPX6

Options (special article):

- Customised SW

It is possible to customise the SW for UBL, OOBD, Switch input and the indicators so that it fits to the customer's system.

Please contact your Sales Engineer at LINAK A/S for further information.

Usage:


- Compatible with all LINAK OpenBus™ Control Boxes, light & full
- Current consumption:
 - *V bus 8 V :* $= 9,3 \text{ mA.}$
 - *V permanent 40 V :* $= 2,3 \text{ mA.}$
- Approvals: LINAK and UL International DEMKO A/S are currently in the process of certifying MJBX0X, re. UL60601 and EN60601-1.

Precautions:

- Always use locking mechanism and O-ring
 - Sockets not used must be fitted with blind plugs to ensure IP Degree
 - When using modular plug with open end, it is up to the customer to maintain the IP degree.
 - The LED end of the UBL cable (0964135) must be mounted on the bed with metal screws, in order to maintain ESD protection.
 - When the MJB with UBL, switch input is being used on a system with battery, the functionality will follow the power down mode of the Control box, see table
- | CB Power mode | | MJB with UBL, OOB & Switch |
|---------------|------------|----------------------------|
| On mains | | Function ok |
| On battery | Power down | No function |
| | "wake up" | Function ok |
- Do not use 2 MJB variants with same device ID on the OpenBus™. This will cause conflicts and the SDT is not able to identify the different products attached.
 - Before the final functional test in the production, it is important that the system is re-powered. This is to make sure that all items have been detected on the OpenBus™.
 - We recommend the end user to make a regular test procedure, in order to prevent failures and hazardous situations on the system, e.g. defective OOB sensor, squeezed cables. The MJB is not able to detect defective 3rd party products.
 - LINAK only takes responsibility for LINAK products, not 3rd party products. Please pay attention to the "Patient Environment" Clause 3.79 - EN60601-1 3rd edition:
Any volume in which contact can occur between a patient and parts of the medical equipment or between a patient and other persons touching parts of the equipment. Justified practice here is 1.5 m. This volume must be subject to the Risk Analysis. It is important to inform the customer about this.
 - There can be a risk of conflict with other accessories such as HB, ACx, SLS etc. it is recommended to make a system overview.

Cable ordering nos.:

Under Bed Light Cable, 0964135-2500, please use metal screws for mounting.

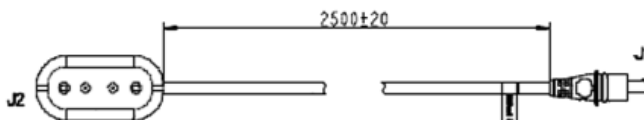

IP Protection degree:

All LINAK products can be cleaned as follow according to their IP protection, which is stated on the product label.

The IP code specifies the degrees of protection provided by the enclosures. Only the protection against ingress of water (second characteristic numeral) is specified, ingress of solid foreign objects (first characteristic numeral) is not specified and therefore replaced by the letter X in the code.

IP Protection	Cleaning instructions
IPX0	Clean with a damp cloth
IPX1	Clean with a damp cloth
IPX2	Clean with a damp cloth
IPX3	Clean with a damp cloth
IPX4	Clean with a damp cloth
IPX5	Wash with a brush and water, but not water under pressure
IPX6	X Wash with a brush and water. The water can be under pressure, but the system must not be hosed down directly with a high pressure cleaner. Max. 20°C
IPX6 Washable	Clean by the use of water tunnels

The materials are resistant to the majority of cleaners and disinfectants used in the hospital and nursing home sector.

Cleaning with a steam cleaner, for example, is not permitted, as it will not be possible to maintain a minimum distance of 300 mm from the electrical parts.

The following directions regarding cleaners and disinfectants must be complied with:

- They must not be highly alkaline or acidic (pH value 6-8).
- They must not contain caustic agents.
- Their contents must not be able to change the structure of the surface or adhesion of the plastic.
- They must not break down grease.

The following disinfectants/detergents are recommended:

- Neodischer Dekonta, dosage 1-3%, made by CHEMISCHE FABRIK Dr. WEIGERT.
- CCOTRADE RW, dosage 0.5%, made by CCOTRADE GmbH.

**Warning!**

The systems must not be sprayed directly with a high-pressure cleaner.

**Warning!**

Interconnecting cables must remain plugged in during cleaning to prevent the ingress of water.

**Precaution!**

When cleaned water must never be sprayed upwards from underneath the bed, because it may enter the cable outlet.


Environmental conditions:

Environmental conditions	
Operating:	
Temperature	5°C to 40°C
Relative humidity	20% to 90% @ 30°C – not condensing
Atmospheric pressure	700 to 1060 hPa
Storage:	
Temperature	-10°C to +50°C
Relative humidity	20% to 90% @ 30°C – not condensing
Atmospheric pressure	700 to 1060 hPa

Valid where nothing otherwise is stated under the specific products in a later section.

Modular plug with wide slit

To be used with **LINAK products**.
Can be connected to all ports in the MJB,
both ports with narrow and wide slit.

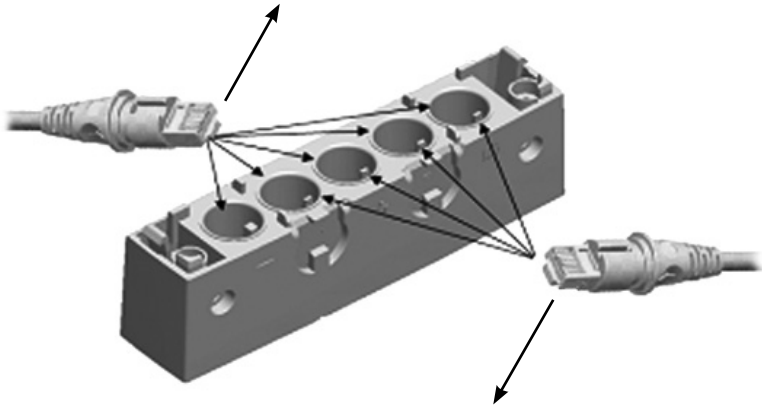

Modular plug with narrow slit

To be used with **3rd party products**.
Can only be connected to ports in the MJB with
narrow slit. This is to prevent 3rd party products
to interfere with the OpenBus™ connections.

Key to symbols

The following symbols are used on the label on the LINAK products.



Type B equipment, as per EN 60601-1

IPXX

Protection against contact/foreign matter (first character) and water (second character) as per EN60529



Class 2 equipment



For indoor use



Double insulated transformer



Protective earth



Alternating Current



Direct current



Attention, consult accompanying documents




Lock function



Release function



Charge indicator

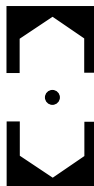

Safety switch/enable button



Demko approval



Fimko approval



APPROVAL
V94265

Australian approval mark

APPROVAL
NO.:97122

Australian approval mark



Recognised - Component Mark



Canadian
Recognised - Component Mark



Recognised
Component Mark for Canada
and the United States



T-Mark



RW-TUV approval



TÜv.Produkt Service



TÜVRheinland



ETL



C-ETL

CS95145V

LGA

LGA



UL Listing Mark



C-UL Listing Mark



C-UL US Listing Mark



UL Listing Mark

LISTED

File E97199

UL file number

File E175209

UL file number

File E151104

UL file number



CSA



PSE-Mark



Product with a thermofuse



For indoor use (House).



Safety isolating transformer.



Electronics scrap



Equipment KL.2 (Double square)



Patient part of type B (Mand)



Patient part of type BF



Earth protection



Equipment class1.

KL.1



Earth



CE Mark




C-TICK

DRAWING APPENDIX

MJB - Modular Junction Box

MJB

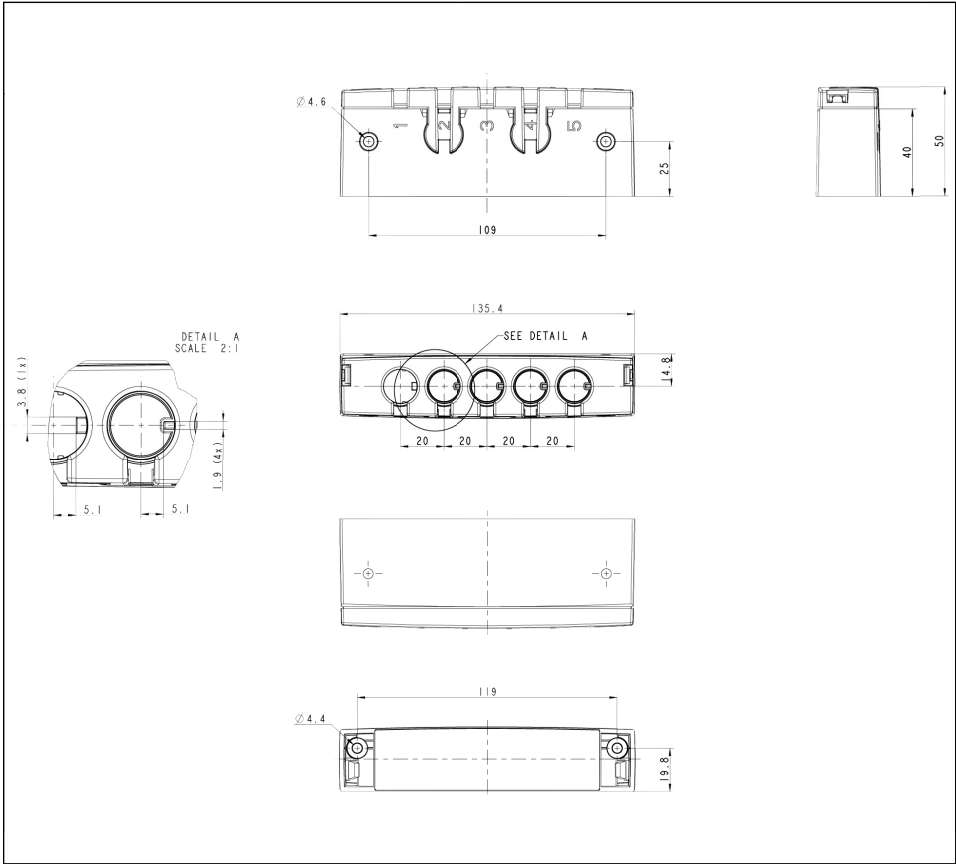

LINAK APPLICATION POLICY

The purpose of the application policy is to define areas of responsibilities in relation to applying a LINAK product defined as hardware, software, technical advice, etc. related to an existing or new customer application.

LINAK products as defined above are applicable for a wide range of applications within the Care and Health, Comfort furniture, Desk and Industry areas. Yet, LINAK cannot know all the conditions under which LINAK products will be installed, used, and operated, as each individual application is unique.

The suitability and functionality of the LINAK product and its performance under varying conditions (application, vibration, load, humidity, temperature, frequency, etc.) can only be verified by testing, and shall ultimately be the responsibility of the LINAK customer using any LINAK product.

It is also the responsibility of the LINAK customer to make and supply a comprehensive user manual of the application.

LINAK shall be responsible solely that the LINAK products comply with the specifications set out by LINAK and it shall be the responsibility of the LINAK customer to ensure that the specific LINAK product can be used for the application in question.

FACTORIES

CHINA

LINAK (Shenzhen) Actuator Systems, Ltd.
Phone: +86 75 58 61 06 656
Fax: +86 75 58 61 06 990
E-mail: sales@linak.cn
www.linak.cn

DENMARK

LINAK A/S - Group Headquarters Guderup
Phone: +45 73 15 15 15
Fax: +45 74 45 80 48
Fax (Sales): +45 73 15 16 13
E-mail: info@linak.com
www.linak.com

SLOVAKIA

LINAK Slovakia s.r.o.
Phone: +421 51 75 63 414
Fax: +421 51 75 63 410
E-mail: jp@linak.sk
www.linak.com

USA

LINAK U.S. Inc.
North and South American Headquarters
Phone: +1 50 22 53 55 95
Fax: +1 50 22 53 55 96
E-mail: info@linak-us.com
www.linak-us.com

SUBSIDIARIES

AUSTRALIA

LINAK Australia Pty. Ltd.
Phone: +61 38 79 69 777
Fax: +61 38 79 69 778
E-mail: sales@linak.au
www.linak.au

AUSTRIA

LINAK Repräsentanz Österreich (Wien)
Phone: +43 (1) 890 74 46
Fax: +43 (1) 92 30 81 09
E-mail: info@linak.de
www.linak.at

BELGIUM & LUXEMBOURG

LINAK Actuator-Systems NV/SA
Phone: +32 (0) 92 30 81 09
Fax: +32 (0) 92 30 88 80
E-mail: info@linak.be
www.linak.be

BRAZIL

LINAK Do Brasil Comércio De Atuadores Ltda.
Phone: +55 (11) 28 32 70 70
Fax: +55 (11) 28 32 70 60
E-mail: info@linak.com.br
www.linak.com.br

CANADA

LINAK Canada Inc.
Phone: +1 50 22 53 55 95
Fax: +1 41 62 55 77 20
E-mail: info@linak.ca
www.linak-us.com

CZECH REPUBLIC

LINAK C&S S.R.O.
Phone: +420 58 17 41 814
Fax: +420 58 17 02 452
E-mail: ponizl@linak.cz
www.linak.cz

DENMARK

LINAK Danmark A/S
Phone: +45 86 80 36 11
Fax: +45 86 82 90 51
E-mail: linak@linak-silkeborg.dk
www.linak.dk

FINLAND

LINAK OY
Phone: +358 10 84 18 700
Fax: +358 10 84 18 729
E-mail: linak@linak.fi
www.linak.fi

FRANCE

LINAK France E.U.R.L.
Phone: +33 (0) 2 41 36 34 34
Fax: +33 (0) 2 41 36 35 00
E-mail: linak@linak.fr
www.linak.fr

GERMANY

LINAK GmbH
Phone: +49 60 43 96 55 0
Fax: +49 60 43 96 55 60
E-mail: info@linak.de
www.linak.de

INDIA

LINAK A/S India Liaison Office
Phone: +91 12 04 39 33 35
Fax: +91 12 04 27 37 08
E-mail: info@linak.in
www.linak.in

IRELAND

LINAK UK Limited - Ireland
Phone: +44 (0) 121 544 2211
Fax: +44 (0) 121 544 2552
E-mail: sales@linak.co.uk
www.linak.co.uk

ITALY

LINAK Italia S.r.l.
Phone: +39 02 48 46 33 66
Fax: +39 02 48 46 82 52
E-mail: info@linak.it
www.linak.it

JAPAN

LINAK K.K.
Phone: +81 45 53 30 802
Fax: +81 45 53 30 803
E-mail: linak@linak.jp
www.linak.jp

MALAYSIA

LINAK Actuators Sdn. Bhd.
Phone: +60 42 10 65 00
Fax: +60 42 26 89 01
E-mail: info@linak-asia.com
www.linak.my

NETHERLANDS

LINAK Actuator-Systems B.V.
Phone: +31 76 5 42 44 40
Fax: +31 76 5 42 61 10
E-mail: info@linak.nl
www.linak.nl

NEW ZEALAND

LINAK New Zealand Ltd.
Phone: +64 95 80 20 71
Fax: +64 95 80 20 72
E-mail: nzsales@linak.com.au
www.linak.co.nz

NORWAY

LINAK Norge AS
Phone: +47 32 82 90 90
Fax: +47 32 82 90 98
E-mail: info@linak.no
www.linak.no

POLAND

LINAK Polska
Phone: +48 (22) 500 28 74
Fax: +48 (22) 500 28 75
E-mail: dkeh@linak.dk
www.linak.pl

REPUBLIC OF KOREA

LINAK Korea Ltd.
Phone: +82 (0) 2 6231 1515
Fax: +82 (0) 2 6231 1516
E-mail: scully@linak.kr
www.linak.kr

RUSSIAN FEDERATION

000 LINAK
Phone: +7 49 52 80 14 26
Fax: +7 49 56 87 14 26
E-mail: info@linak.ru
www.linak.ru

SPAIN

LINAK Actuadores, S.L.u
Phone: +34 93 58 82 777
Fax: +34 93 58 82 785
E-mail: linakact@linak.es
www.linak.es

SWEDEN

LINAK Scandinavia AB
Phone: +46 87 32 20 00
Fax: +46 87 32 20 50
E-mail: info@linak.se
www.linak.se

SWITZERLAND

LINAK AG
Phone: +41 43 38 83 188
Fax: +41 43 38 83 187
E-mail: info@linak.ch
www.linak.ch

DISTRIBUTORS

ARGENTINA

Novotec Argentina SRL
Phone: +54 (11) 4303 8900 / 89
Fax: +54 (11) 4032 0184
E-mail: info@novotecargentina.com
www.novotecargentina.com

AUSTRALIA

Ballarat Industrial Supplies
www.ballind.com.au
BL Shipways & Co
www.blshipway.com.au
Gas Strut Marine and Industrial
www.gasstrutmarine.com.au
Prime Motion & Control
www.primemhd.com.au
West Vic Industrial Supplies
www.westvicindustrial.com.au

COLOMBIA

MEM Ltda
Phone: +57 (1) 334 7666
Fax: +57 (1) 282 1684
E-mail: servicioalcliente@memltda.com.co
www.memltda.com.co

INDONESIA

Pt. Himalaya Everest Jaya
Phone: +62 21 54 48 956 / 65
Fax: +62 21 61 94 658 / 1925
E-mail: hejplastic-div@centrin.net.id
www.hej.co.id

IRAN

Bod Inc.
Phone: +98 21 88 99 86 35 - 6
Fax: +98 21 88 95 44 81
E-mail: info@bod.ir
www.bod.ir

MEXICO

ILSA S.A. de C.V.
Phone: +52 (55) 53 88 39 60
Fax: +52 (55) 53 88 39 66
E-mail: linak@ilsamexico.com
www.ilsamexico.com

PERU

Percy Martín Del Aguila Ubillus
Phone: +51 99 88 39 879

RUSSIAN FEDERATION

000 FAM
Phone: +7 81 23 31 93 33
Fax: +7 81 23 27 14 54
E-mail: purchase@fam-drive.ru
www.fam-drive.ru

SINGAPORE

Servo Dynamics Pte. Ltd.
Phone: +65 68 44 02 88
Fax: +65 68 44 00 70
E-mail: servodynamics@servo.com.sg
www.servo.com.sg

SOUTH AFRICA

Industrial Specialised Applications CC
Phone: +27 11 31 22 292 or
+27 11 20 77 600
Fax: +27 11 31 56 999
E-mail: garth@isagroup.co.za
www.saza.co.za

UNITED ARAB EMIRATES

Mechatronics
Phone: +971 4 26 74 311
Fax: +971 4 26 74 312
E-mail: mechtron@emirates.net.ae
www.mechtronics.ae

For contact details on other countries please visit www.linak.com or contact:

LINAK INTERNATIONAL

Phone: +45 73 15 15 15
Fax: +45 74 45 90 10
E-mail: info@linak.com
www.linak.com

Terms of use

The user is responsible for determining the suitability of LINAK products for specific application. LINAK takes great care in providing accurate and up-to-date information on its products. However, due to continuous development in order to improve its products, LINAK products are subject to frequent modifications and changes without prior notice. Therefore, LINAK cannot guarantee the correct and actual status of said information on its products. While LINAK uses its best efforts to fulfil orders, LINAK cannot, for the same reasons as mentioned above, guarantee the availability of any particular product. Therefore, LINAK reserves the right to discontinue the sale of any product displayed on its website or listed in its catalogues or other written material drawn up by LINAK. All sales are subject to the Standard Terms of Sale and Delivery for LINAK. For a copy hereof, please contact LINAK.



WE IMPROVE YOUR LIFE